

This Page Is Inserted by IFW Operations  
and is not a part of the Official Record

## **BEST AVAILABLE IMAGES**


Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning documents *will not* correct images,  
please do not report the images to the  
Image Problem Mailbox.**

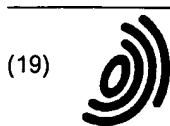

Europäisches Patentamt  
European Patent Office  
Office européen des brevets



(11) **EP 0 992 795 A1**

(12)

## EUROPEAN PATENT APPLICATION

(43) Date of publication:  
12.04.2000 Bulletin 2000/15

(51) Int Cl.7: **G01N 33/70**

(21) Application number: **99117955.7**

(22) Date of filing: **15.09.1999**

(84) Designated Contracting States:  
**AT BE CH CY DE DK ES FI FR GB GR IE IT LI LU  
MC NL PT SE**  
Designated Extension States:  
**AL LT LV MK RO SI**

(30) Priority: **28.09.1998 US 161768**

(71) Applicant: **Bayer Corporation**  
**Pittsburgh, PA 15205-9741 (US)**

(72) Inventors:  
• **Cast, Todd K.**  
**Granger, Indiana 46530 (US)**  
• **Pugia, Michael J.**  
**Granger, Indiana 46530 (US)**

(74) Representative: **Burkert, Frank et al**  
**Bayer AG,**  
**Konzernbereich RP**  
**Patente und Lizenzen**  
**51368 Leverkusen (DE)**

(54) **Improved method for the detection of creatinine**

(57) Disclosed is an improvement to the assay for creatinine in fluid test samples by contacting the test sample with a reagent formulation containing cupric ions, a hydroperoxide and an oxidizable dye. The improvement involves the inclusion of one or more selected quinolines in the reagent formulation. When the qui-

noline is substituted with methyl in the 2-5 position and/or a methoxy in the 2-6 position autooxidation of the oxidizable dye is reduced. When the quinoline is further substituted in an otherwise unsubstituted carbon atom in the 2-8 position, the assay's resistance to interference from hemoglobin and ascorbate is increased.

EP 0 992 795 A1

**Description****Background of the Invention**

5 [0001] The present invention is an improvement to the method of detecting creatinine in samples of body fluid. There is disclosed in U.S. Patent 5,374,561 a method for the detection of creatinine in urine which involves contacting a urine sample suspected of containing creatinine with cupric ions, a hydroperoxide and an oxidizable dye which provides a colored response in the presence of oxygen free radicals and a hydroperoxide. In this assay, the first step involves the formation of a  $\text{Cu}^{++}$ -creatinine chelated complex. The oxidizable dye is oxidized by the transfer of an electron therefrom to the  $\text{Cu}^{++}$ -creatinine complex to provide the non-reactive  $\text{Cu}^{+}$ -creatinine form which is regenerated to  $\text{Cu}^{++}$ -creatinine by the loss of an electron to the hydroperoxide. This assay works well in the absence of hemoglobin and ascorbate, but in the presence of hemoglobin and/or ascorbate, which are normally present in clinical urine samples, the precision of the assay is affected due to the tendency of hemoglobin and/or ascorbate to oxidize the dye which results in false positive results. In addition, the redox dye can auto-oxidize over time in the presence of copper ion thereby reducing the shelf life of this type of assay. Accordingly, it would be desirable and it is an object of the present invention to provide an improved reagent system for the detection of creatinine in a fluid test sample in which the tendency of hemoglobin and/or ascorbate in the test sample to cause falsely positive results and the auto-oxidation of the redox dye is reduced.

**Summary of the Invention**

20 [0002] The present invention is an improvement to the method for the detection or determination of creatinine in a fluid test sample which involves contacting the test sample with cupric ions, a hydroperoxide and an oxidizable dye which provides a detectable colored response upon being oxidized when the cupric ions and creatinine combine to form a complex having peroxidase activity. This type of assay suffers from two drawbacks. First of all, the oxidizable dye can undergo auto-oxidation in the presence of transition state metals thereby shortening the assay reagents' shelf life. In addition, the presence of hemoglobin and/or ascorbate in the fluid test sample can cause false positive results due to their ability to oxidize the oxidizable dye in the absence of creatinine. It has now been discovered that these problems can be ameliorated by the addition of certain quinoline compounds to the assay reagents. Typically, those quinolines which have a methyl group in the 2-5 position and/or methoxy in the 2-6 position have been found to reduce or eliminate the problem of auto-oxidation. Certain quinoline compounds, i.e. those which have at least one nitro, amine, halogen or hydroxyl group in the 2-8 position and a methyl or methoxy group in position 2-6, have been found to have a positive effect on the problem of auto-oxidation and also to substantially reduce the tendency of hemoglobin and/or ascorbate to oxidize the dye thereby reducing the incidence of false positive results in the assay.

**Description of the Invention**

35 [0003] A novel method for the colorimetric determination of creatinine in fluid test samples is disclosed in U.S. Patent 5,374,561. This method takes advantage of the peroxidase activity of creatinine/cupric ion metal complexes. In this assay, the chelation of cupric ions by creatinine produces a complex having peroxidase like activity which is capable of catalyzing the oxidation of an oxidizable dye such as tetramethyl benzidine by a hydroperoxide such as diisopropyl benzene dihydroperoxide to thereby provide a colored response. It has been discovered more recently that the accuracy of this assay is adversely affected by auto-oxidation of the oxidizable dye in the presence of metals which catalyze oxidation of redox indicators by an oxygen source through electron transfer, thereby reducing the reagent's shelf life. Of greater concern is the presence of hemoglobin and/or ascorbate in the test fluid which can oxidize the dye in the absence of creatinine resulting in false positive results.

40 [0004] Traylor et al discuss the capability of aromatic pyridines to increase peroxidase activity in Tetrahydron, 40: 1984; 553-68. This teaching suggests that aromatic amines like quinoline cause redox dye oxidation in the presence of hemoglobin. In U.S. Patent 3,853,472 to Rittersdorf there is disclosed a method for the detection of blood in urine which indicates that quinolines increase peroxidase activity as evidenced by increased redox dye oxidation.

50 [0005] It has now been discovered that certain quinoline compounds, rather than increasing peroxidase activity, can reduce the problems associated with auto-oxidation of the oxidizable dye. Those quinolines which are useful for this purpose are characterized by being substituted with methyl in the 2-5 position and/or methoxy in the 2-6 position. Thus, by adding one or more of these quinoline compounds to the reagent formulation at a concentration of from 10 to 300 mM, preferably from 50 to 150 mM, the shelf life of the reagent system can be extended.

55 [0006] Another, more limited class of quinolines can be added to the  $\text{Cu}^{++}$ -creatinine reagent system to control the adverse effects of hemoglobin and/or ascorbate on the oxidizable dye. These quinolines, which also have the effect of reducing the effects of auto-oxidation are characterized by at least one nitro, amine, halogen or hydroxyl group in the 2-8 position of the quinoline and a methyl in position 2-5 or methoxy group in position 2-6. They are typically added

to the assay formulation in amounts as described above.

[0007] In the copper<sup>++</sup> assay the source of cupric ion may be any soluble copper salt whose anion does not detrimentally interact with the reaction for the colorimetric detection of creatinine. Suitable salts include copper sulfate, nitrate oxide, hydroxide, phosphate, iodide, chloride, bromide, acetate or oxalate. Other soluble cupric salts may be used provided that they allow formation of the Cull-Creatinine complex. Those salts whose anion binds too strongly to the copper will not allow the copperII-Creatinine complex to be formed. Accordingly, Cull complexes such as those formed between cupric ions and EDTA, HEDTA, EGTA and DTPA would not release sufficient Cull for formation of the Cull-Creatinine complex. It has been observed that the citrate and sulfate salts have the lowest blank reactivity and, accordingly, they are preferred. Cupric citrate is particularly preferred due to its exhibiting the least blank reactivity and the greatest formation of the Cull-Creatinine complex. Salts which oxidize the dye in the absence of creatinine are less desirable. Salts such as cupric 2,2'-bipyridine can cause a significant oxidation of TMB in the absence of creatinine, and are, therefore, unsuitable for use in the present invention. When copper citrate is used as the cupric ion source, the concentration of citrate ion should be at least that of copper. An excess of citrate ion of at least twice that of the copper ion is preferred to ensure complete complexation of Cull by the citrate.

[0008] Typically, when urine is the aqueous fluid being tested, the concentration of cupric ion will be from 5 to 30 mM since the reference range of creatinine in urine is 3 to 20 mM. This range would vary in other fluids such as serum where one would preferably employ a concentration of cupric ion in the range of from 0.05 to 0.30 mM. The Cuprous ion tends to cause some background interference due to oxidation of the dye in the absence of creatinine. Accordingly, Cu<sup>+</sup> salts cannot be used.

[0009] Suitable oxidizable indicators include, for example, benzidine; o-tolidine; a 3,3',5,5'-tetraalkylbenzidine wherein the alkyl group includes from one to about six carbon atoms; o-dianisidine; 2,7-diaminofluorene; bis-(n-ethylquinol-2-one)-azine; (N-methylbenzthiazol-2-one)-(1-ethyl-3-phenyl-5-methyltriazol-2-one)-azine or combinations thereof.

[0010] Suitable hydroperoxides for use in the present invention include cumene hydroperoxide; 5-butyl hydroperoxide; diisopropylbenzene hydroperoxide; 1-hydroxycyclohexane-1-hydroperoxide; 2,5-dimethyl-hexane-2,5-dihydroperoxide; paramenthane hydroperoxide; 1,4-diisopropylbenzene monohydroperoxide; p-t-butyl-isopropylbenzene hydroperoxide; 2-( $\alpha$ -hydroperoxyiso-propyl)-6-isopropyl-naphthalene; tetralin hydroperoxide or combinations thereof.

[0011] Typically, the reagent system, comprising the soluble copper salt, hydroperoxide and oxidizable indicator will be dissolved in water. However, organic solvents can be incorporated into the system provided they do not interfere with the assay mechanism. The concentration of the hydroperoxide and oxidizable indicator will normally range from 10 to 150 mM with a range of from 60 to 100 mM being preferred.

[0012] In the practice of the invention, the assay can be performed in either the wet or the dry (test strip) format. In carrying out the assay, the test sample is mixed with the copper salt, e.g. cupric citrate, the dye and the hydroperoxide at a buffered pH, preferably from 4.0 to 9.0, through the use of a reagent strip or aqueous and acetonitrile solutions of reagents. Reagent strips are prepared in the conventional manner of dipping an absorbant carrier into an aqueous solution of the cupric salt and buffers, drying the carrier and then dipping it into an organic solution of the dye and hydroperoxide with subsequent drying.

[0013] The present invention is further illustrated by the following Examples:

#### Example I

[0014] A creatinine reagent for use in this experiment was prepared by a two dip method in which a strip (1/5 inch x 1/5 inch) of Whatman cellulose filter paper grade 3 MM pad on a 1/5 inch by 1/2 inch polystyrene handle was dipped into a first dip solution with subsequent drying at 110°C for 10 minutes followed by dipping it into a second dip solution and drying. The dip solutions were formulated as follows:

1st Dip	
Component	Concentration
CuSO <sub>4</sub>	30 mM
Citrate	50 mM
Buffer: Glycerol-2-phosphate	750 mM
Surfactant: Aerosol IB-45	1.5%
pH	6.80
Water	

2nd Dip	
Component	Concentration
3,3',5,5'-Tetramethylbenzidine (TMB)	33 mM
Diisopropyl Benzene Dihydroperoxide (DBDH)	80 mM
Polymer: Plasdane	2.5%
Quinoline; 4-Hydroxy-2-methylquinoline	100 mM
Ethyl Orange Alcohol	0.32%

[0015] This creatinine reagent demonstrates a large stability shift as shown in Table 1.

Table 1

The Stability Shift of the Creatinine Formula	
Storage Temperature	ND result for a 30 mg/dl sample observed within 3 months of storage
25° C	21
30° C	23
40° C	24
50° C	19
60° C	18

[0016] A specimen containing either 30, 100 or 200 mg/dL of creatinine should read the expected result throughout the shelf life of the strip. However, in the stability shift demonstrated by Table 1, numerical designations as high as 24 were observed with a sample containing 30 mg/dL creatinine. This result is unacceptably high since visual results are measured by interpolation using a color chart with numerical designations of 10, 20 and 30 for the colors corresponding to 30, 100 and 200 mg/dL of creatinine concentration. For example, a result for a sample with 250 mg/dL should produce a numerical designation result of 25 or halfway between the results expected for samples containing 100 mg/dL and 200 mg/dL creatinine. As can be determined from Table 1, numerical designations (ND) as high as 24 were observed with a sample containing 30 mg/dL creatinine that should read as 10.

[0017] Several theories were pursued concerning the cause of the stability shift. By elimination of components, it was determined that the problem was related to auto-oxidation of the redox indicator, and it was discovered that when lepidine (4-methylquinoline) 100 mM was added to the creatinine formula, the visual stability of the assay improved dramatically as demonstrated by the data of Table 2.

Table 2

The Stability Shift of the Creatinine Formula with Lepidine	
Storage Temperature	ND result for a 30 mg/dl sample observed within 3 months of storage
25° C	11
30° C	12
40° C	10
50° C	11
60° C	10

[0018] As can be determined from Table 2, numerical designations (ND) of 10 to 12 were observed with a sample containing 30 mg/dL which would ideally read as 10.

[0019] While better stability was achieved with the lepidine containing formulation, there were also presented several problems. The addition of lepidine decreased the dynamic range of the test from 30 to 300 mg/dL to about 30 to 200 mg/dL. The presence of lepidine caused the creatinine formulation to turn into an occult blood test since the purpose of the lepidine is to enhance the sensitivity of the chemistry to hemoglobin. The presence of 1 mg/dL hemoglobin caused the lepidine containing formulation to read as a formulation containing 30 mg/dL creatinine as more than 300 mg/dL. This amount of interference due to hemoglobin is unacceptable for use in clinical urines which often contain blood.

[0020] The screening of other quinoline compounds similar to lepidine was carried out in order to find a derivative solving the stability problem without causing hemoglobin sensitivity and reduced dynamic range. Thirty quinoline de-

rivatives were screened at concentrations of 25 and 100 mM. A list of the quinolines is set out in Table 3.

Table 3

Quinoline Derivatives Screened in Creatinine Formulation	
Quinoline Compounds Tested	
8-Ethoxyquinoline sulfonic acid	8-Nitroquinoline
2-Chloroquinoline	6-Methylquinoline
2-Chlorolepidine	7-Methylquinoline
5-Chloro-8-hydroxyquinoline	3-Aminoquinoline
2-Hydroxyquinoline	
8-Hydroxyquinoline	
6-Methoxyquinoline	
3-Methylquinoline	
8-Methylquinoline	4-Hydroxy-2-methylquinoline
5-Nitroquinoline	
6-Nitroquinoline	3-Bromoquinoline
8-Nitroquinoline	
Quinaldine (2-Methylquinoline)	
4-Chloro-quinaldine	2-quinolinecarbonitrile
8-Hydroxy-quinaldine	3-quinolinecarbonitrile

[0021] The quinolines selected for evaluation were similar to lepidine (4-methylquinoline) but with the methyl group at another position on the quinoline ring, had an amine, hydroxy or nitro group on the quinoline ring or were substituted at positions 2 or 8 since it was hypothesized that substitution at the 2 or 8 position would be effective for reduction of hemoglobin interference due to steric factors. Each quinoline derivative was incorporated into the formulation described in Example I at a concentration of 25 and 100 mM and examined for formulation compatibility (dip compatibility, impregnated paper appearance and color development); tested for hemoglobin sensitivity at 1 mg/dL hemoglobin and evaluated for its effect on thermal stability.

[0022] Of the quinolines tested, twelve were found to provide hemoglobin resistance which was better than the creatinine formula without quinoline. These quinolines were 2-Chlorolepidine, 2-Chloroquinoline, 2-Hydroxyquinoline, 8-Methylquinoline, 8-Nitroquinoline, 3-Aminoquinoline, 4-Chloroquinoline, 5-Nitroquinoline, 6-Nitroquinoline, 4-Hydroxy-2-methylquinoline, 8-Hydroxyquinoline and 3-Bromoquinoline. Comparison of these twelve quinolines structurally supports the hypothesis that the substitution of a halogen, nitro, hydroxyl or amine group on the ring in the 2 or 8 position is the most effective means of preventing hemoglobin interference with the creatinine reagent.

[0023] Only three quinolines; 6-Methoxyquinoline, 3-Methylquinoline and 4-Hydroxy-2-methylquinoline demonstrated a stability improvement in the creatinine formulation similar to that of lepidine suggesting that a methyl in the 2, 3, 4 or 5 position or methoxy in the 2, 3, 4, 5 or 6 position is required to provide the requisite stability improvement. Two of the three quinolines, 6-Methoxyquinoline and 3-Methylquinoline gave stability and reactivity performance similar to that of lepidine. These compounds increased the overall reactivity and reduced the dynamic range of the test. The other compound, 4-Hydroxy-2-methylquinoline, did not exhibit performance similar to lepidine since it gave reactivity similar to the control with even greater resolution between levels. It was also the only quinoline compound that was found to provide hemoglobin resistance and the desired stability improvement. This quinoline is less hazardous than lepidine and is odorless. Based on the significant improvements in stability and hemoglobin resistance, 4-Hydroxy-2-methylquinoline is the preferred quinoline for use in the present invention. Quinolines substituted with amino, nitro, hydroxyl or halogen in the 2-8 position, methoxy in the 2-6 position and/or methyl in the 2-5 position would provide similar improvement due to their ability to interfere with the quinoline's ability to coordinate with hemoglobin.

[0024] The concentration of 4-Hydroxy-2-methylquinoline was varied from 50 mM to 120 mM in 10 mM increments to determine if a lower level would provide better high temperature stability. A definite relationship was observed between quinoline concentration and thermal stability with the apparent optimal quinoline concentration for this formulation being 80-110 mM. As can be determined from Table 4, numerical designations (ND) observed with samples containing 10 to 300 mg/dL creatinine were closer to the expected ND with the improved formula. The ND values observed with the quinoline additive were not affected by storage at 40°C for 4 weeks while the ND values observed with the original formulations were affected. The analytical range of the improved formulation was also greater.

Table 4

Creatinine Improvement Formulation versus Original Creatinine Reagent												
Concentration			10 mg/dl	50 mg/dl	100 mg/dl	200 mg/dl	300 mg/dl					
Expected ND			10	20	30	40	50					
Formula	Storage temp.	Weeks stored	observed mean	ND sd	observed mean	ND sd	observed mean	ND sd	observed mean	ND sd	observed mean	
Original Creatinine with no quinoline	-20C	2	7.8	1.4	9.8	1.9	20.3	1.3	29.8	1.6	40.0	2.8
		4	7.3	0.5	12.2	1.6	23.7	2.4	32.7	2.3	43.0	2.3
	40C	2	9.2	1.4	14.3	2.5	23.2	1.9	30.1	2.0	40.6	1.9
		4	12.0	1.1	18.3	2.3	24.9	2.6	32.5	3.2	42.8	3.7
Improved Creatinine with 4-hydroxy 2-methylquinoline	-20C	2	8.1	1.6	17.3	2.0	27.2	2.5	37.9	2.7	45.5	2.3
		4	8.5	1.6	18.7	1.8	26.3	2.5	40.6	3.3	45.3	2.1
	40C	2	10.4	0.7	20.5	2.4	28.1	2.4	38.4	3.2	46.3	2.0
		4	11.3	1.2	19.0	4.0	25.8	4.1	36.5	4.0	44.6	4.0

[0025] Interference testing was carried out with ascorbic acid and hemoglobin. Concentrations of 1 to 10 mg/dL hemoglobin and 25 to 440 mg/dL ascorbic acid were tested in a 50 mg/dL creatinine urine standard with both the improved creatinine formulation and the original creatinine reagent. Results of interference testing, using a CLINITEK®-50 reflectance spectrometer, are set out in Table 5.

Table 5

Hemoglobin and Ascorbate Resistance of Current and Improved Creatinine Reagents									
Form.	[Interferent] (mg/dL)	Stab Cond.	CLINITEK® 50 Decode Result						
			<10mg/dL		50mg/dL			100mg/dL	
			Mean	SD	Mean	SD	% Change	Mean	SD
Original creatinine formula without quinoline	none	2W25C	541	12	486	10		383	7
	1 Hb	2W25C			489	10	-0.3		
	2 Hb	2W25C			461	11	24.3		
	4Hb	2W25C			415	8	68.9		
	6 Hb	2W25C			382	7	101.0		
	10 Hb	2W25C			335	12	146.6		
	25 AA	2W25C			395	12	88.4		
Improved creatinine	440 AA	2W25C			284	12	196.1		
	none	2W25C	528	9	437	7		327	9
	1 Hb	2W25C			443	13	-0.5		
	2 Hb	2W25C			427	6	9.1		
	4Hb	2W25C			418	9	17.3		
	6 Hb	2W25C			402	10	31.8		
	10 Hb	2W25C			380	7	51.8		
	25 AA	2W25C			429	6	7.3		
	440 AA	2W25C			442	12	-0.5		

1) Hb = hemoglobin

2) AA = ascorbate

3) Decode = a number representing the reflectance of color off the reagent as measured by the CLINITEK® 50 instrument. The lower the number, the more color has been generated. The typical decode responses for three creatinine levels; 10, 50 and 100 mg/dL; are provided for comparison.

[0026] Referring to Table 5, the standard creatinine formulation was found to exhibit both greater hemoglobin and ascorbate sensitivity than the improved formulation of the present invention. The test, absent the quinoline of the present invention, shows a 147% increase in reactivity with 10 mg/dL hemoglobin concentration compared to 52% for the quinoline containing formulation. In general, a three fold improvement in hemoglobin resistance was observed with the improved creatinine test at all concentrations of hemoglobin tested. In actual practice, the urine sample being tested would not contain such a high concentration of hemoglobin because the devices used for this sort of assay are labeled




with a direction not to use them on visibly bloody urine which occurs at about 5 mg/dL hemoglobin.

[0027] With regard to ascorbate resistance, the current creatinine test exhibited an 82% gain in reactivity in the presence of 25 mg/dL ascorbate. In comparison, the improved creatinine test demonstrated resistance to 25 mg/dL ascorbate and also to ascorbate at a concentration of 440 mg/dL, twenty times the amount normally found in clinical urines. Thus, both hemoglobin and ascorbate resistance are significantly improved with the reagent of the present invention.

## Claims

1. An assay for creatinine in a fluid test sample which assay comprises contacting the fluid test sample suspected of containing creatinine with cupric ions, a hydroperoxide and an oxidizable dye, characterized in the assay comprises inclusion in the fluid test sample of a substituted quinoline compound which inhibits auto-oxidation of the oxidizable dye in the presence of cupric ion.

2. The assay of Claim 1 wherein the substituted quinoline is of the formula:

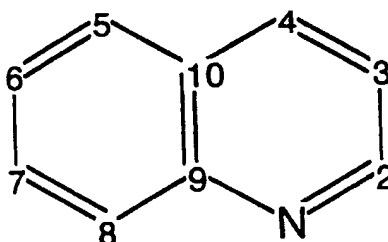

which is substituted with a methyl group at a carbon atom occupying position 2-5 of the molecule and/or with a methoxy group at an otherwise unsubstituted carbon atom occupying position 2-6 of the molecule.

3. The assay of any of Claims 1 and 2 wherein the quinoline compound is 4-Hydroxy-2-methyl quinoline, 6-Methoxyquinoline or 3-Methylquinoline.

4. The assay of any of Claims 1 to 3 wherein the substituted quinoline is present in the reagent formulation at a concentration of from 10 to 300 mM.

5. The assay of any of Claims 1-4 wherein the concentration is from 50 to 150 mM.

6. The assay of any of Claims 1 to 5 wherein the fluid test sample comprises urine.

7. The assay any of Claims 1 to 6 wherein the reagent formulation is dried onto a strip of absorbant material to provide a dry assay in which the strip is contacted with the fluid test sample.

8. An assay for creatinine in urine in which the urine is contacted with a reagent system comprising cupric ions, a hydroperoxide and an oxidizable dye together with 4-hydroxy-2-methylquinoline.

9. The assay of Claim 8 wherein the 4-hydroxy-2-methylquinoline is present in the reagent formulation at a concentration of from 10 to 300 mM, the hydroperoxide is diisopropyl benzene and the oxidizable dye is 3,3',5,5'-tetramethylbenzidine.

10. The assay of any of Claims 8 and 9 wherein the reagent formulation is dried onto a strip of absorbant material to provide a dry assay in which the strip is contacted with the urine.



European Patent  
Office

# EUROPEAN SEARCH REPORT

Application Number  
EP 99 11 7955

DOCUMENTS CONSIDERED TO BE RELEVANT			
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int.Cl.7)
A	EP 0 814 335 A (BAYER AG) 29 December 1997 (1997-12-29) * the whole document *	1-10	G01N33/70
A	US 4 755 472 A (ISMAIL IBRAHIM A ET AL) 5 July 1988 (1988-07-05) * column 3, line 55 - line 69 *	1-10	
A	US 4 886 760 A (BEN-MICHAEL ABRAHAM) 12 December 1989 (1989-12-12) * column 5, line 27 - line 49 * * column 7, line 9 - line 17 *	1-10	
A	EP 0 043 469 A (MILES LAB) 13 January 1982 (1982-01-13) * the whole document *	1-10	
A	US 5 081 040 A (PATEL CHANDRAVADAN ET AL) 14 January 1992 (1992-01-14) * column 5, line 3 - line 6 * * column 7, line 22 - line 24 *	1-10	
			TECHNICAL FIELDS SEARCHED (Int.Cl.7)
			G01N C12Q
The present search report has been drawn up for all claims			
Place of search THE HAGUE		Date of completion of the search 14 January 2000	Examiner Hart-Davis, J
<p>CATEGORY OF CITED DOCUMENTS</p> <p>X : particularly relevant if taken alone Y : particularly relevant if combined with another document of the same category A : technological background O : non-written disclosure P : intermediate document</p> <p>T : theory or principle underlying the invention E : earlier patent document, but published on, or after the filing date D : document cited in the application L : document cited for other reasons</p> <p>&amp; : member of the same patent family, corresponding document</p>			

EPO FORM 1503 03/82 (P04/001)

**ANNEX TO THE EUROPEAN SEARCH REPORT  
ON EUROPEAN PATENT APPLICATION NO.**

EP 99 11 7955

This annex lists the patent family members relating to the patent documents cited in the above-mentioned European search report. The members are as contained in the European Patent Office EDP file on  
The European Patent Office is in no way liable for these particulars which are merely given for the purpose of information.

14-01-2000

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
EP 0814335 A	29-12-1997	US 5733787 A	31-03-1998
		AU 2496197 A	08-01-1998
		CA 2198957 A	17-12-1997
		JP 10062412 A	06-03-1998
US 4755472 A	05-07-1988	AT 67789 T	15-10-1991
		AU 569436 B	28-01-1988
		AU 6501486 A	23-07-1987
		CA 1290661 A	15-10-1991
		DE 3773218 A	31-10-1991
		DK 20187 A	17-07-1987
		EP 0230229 A	29-07-1987
		FI 870136 A,B,	17-07-1987
		IE 60127 B	01-06-1994
		JP 1772500 C	14-07-1993
		JP 4066314 B	22-10-1992
		JP 62169053 A	25-07-1987
		NO 870020 A,B,	17-07-1987
US 4886760 A	12-12-1989	IL 74204 A	23-12-1990
		CA 1329475 A	17-05-1994
		EP 0196743 A	08-10-1986
		FI 860444 A	01-08-1986
		JP 61182558 A	15-08-1986
EP 0043469 A	13-01-1982	CA 1143634 A	29-03-1983
		JP 57013358 A	23-01-1982
		JP 63067139 B	23-12-1988
US 5081040 A	14-01-1992	NONE	

EPO FORM P0459

For more details about this annex : see Official Journal of the European Patent Office, No. 12/82